

AARON® 2250™

Bovie®

ELECTROSURGERY

FOR THE HOSPITAL OPERATING ROOM AND SURGERY CENTER

DISTRIBUTED BY




SOUTHEASTERN
BIOMEDICAL ASSOCIATES INC

3546 Hickory Blvd. Hudson, NC 28638

888.310.7322 www.sebiomedical.com

Email: info@sebiomedical.com

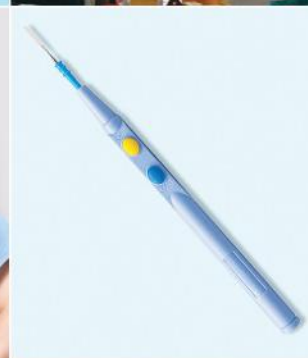

ELECTROSURGICAL GENERATOR

CUTTING EDGE

200 watts of Cutting power

COMPETENCE

Enough power to perform any procedure in a surgery center and most Hospital OR procedures

VERSATILITY

Offers Pure Cut, ten (10) modes of Blend, two (2) modes of Coagulation, Bipolar Coagulation and BovieFCFS™ technology allowing the Doctor to move almost effortlessly from Monopolar to Bipolar during procedure

AFFORDABLE

The highest in safety and quality were built into a 200 watt generator for maximum reliability with economy

SECURITY

The Aaron 2250 is designed with state-of-the-art safety features including Bovie NEM™ pad sensing technology which monitors the return electrode for optimum patient protection. BovieDED™ monitors all functions of the generator to ensure constant safety for the patient, doctor and OR staff





RELIABILITY

The Aaron 2250 is backed by Bovie's two (2) year manufacturer's limited warranty

OTHER OPTIONAL ACCESSORIES

Disposable & Reusable Pencils*	Disposable Grounding Pads*	Monopolar & Bipolar Footswitches
		
ESP1 & ESPR	ESRE & ESREC	BV-1253B & BV-1254B

Disposable Suction Coagulator*	Reusable SuperCut™ Tungsten Needle*	Resistick II™ Coated Blade*
		
SCH10 Handswitch	ES60R	ES01T

Disposable and Reusable Electrodes*		Bipolar Forceps*	Bipolar Forcep Cord*
Blade	Ball	Bayonet	Reusable
			
ES01	ES20	A821	A827V

*More styles available

Industry-standard accessories fit most electrosurgical generators

Additional accessories available — contact your Bovie representative for details.

TECHNICAL SPECIFICATIONS

Technical Specifications...

Dimensions:		Line Voltage:	Line Frequency:
Height	15.30 cm (6.00 in.)	100-240 VAC	50 - 60 Hz
Width	31.10 cm (12.30 in.)		
Depth	41.30 cm (16.30 in.)	Output Type	Duty Cycle
Weight	< 8.8 kg (< 19 lbs)	RF ISOLATED <input checked="" type="checkbox"/> F	10 S/30 S
Monopolar	Output Power	Output Frequency	Rep Rate
Cut	200 W @ 300 Ω	490 kHz ± 5 kHz	N/A
Blend (MAX)	200 W @ 300 Ω	490 kHz ± 5 kHz	30 kHz ± 5 kHz
Coagulation	120 W @ 500 Ω	490 kHz ± 5 kHz	30 kHz ± 5 kHz
Fulguration	80 W @ 500 Ω	490 kHz ± 5 kHz	30 kHz ± 5 kHz
Bipolar	80 W @ 150 Ω	490 kHz ± 5 kHz	30 kHz ± 5 kHz
Warranty: limited two (2) year on unit only			

Aaron® 2250™ Unit



BV-IDS-CS Stand
(Optional)

Smoke Evacuation System
(Optional)



Bovie Medical Corporation • 5115 Ulmerton Road • Clearwater, FL 33760 • www.boviemed.com • sales@boviemed.com
U.S. Phone 1-800-537-2790 • Fax 1-800-323-1640 • Int'l Phone +1-727-384-2323 • Fax +1-727-347-9144

Bovie® and Aaron® are registered trademarks of Bovie Medical Corporation.